


VT300

TILTING WALL STAND

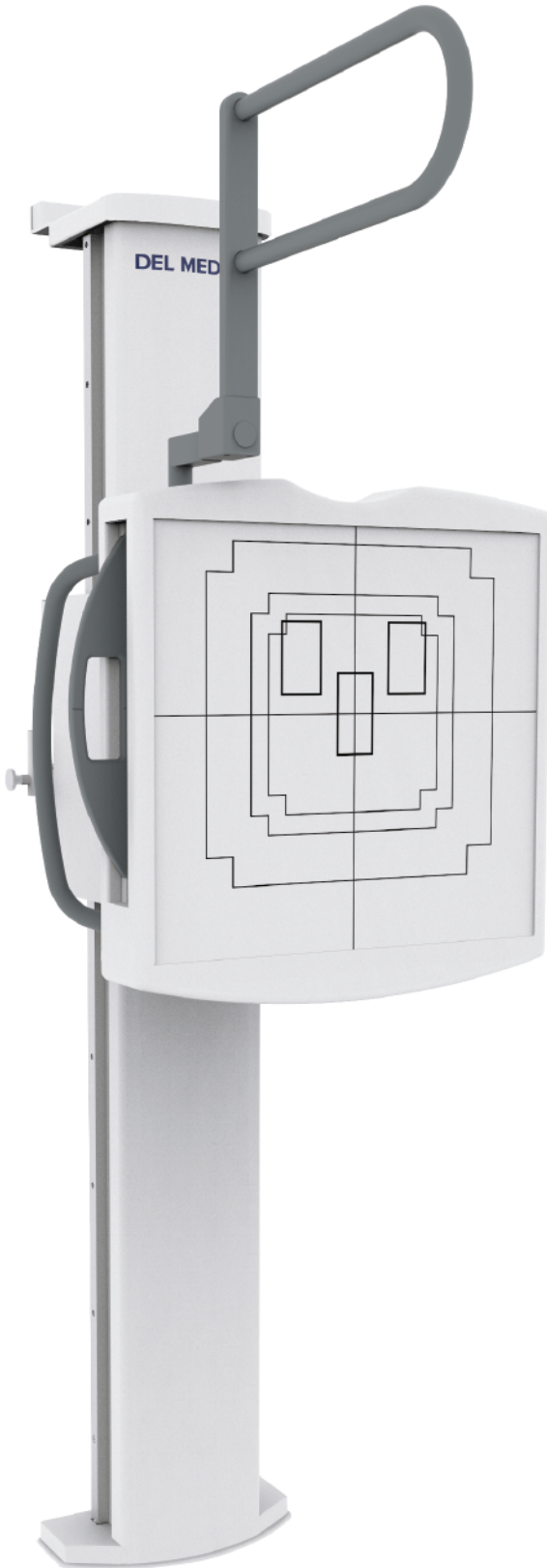

DEL MEDICAL's tilting VT300 wallstand is an easily positioned and secure housing for imaging receptors such as film, CR and all cassette sized digital flat panel detectors. The image receptor housing has a tilt range of $+90^{\circ}/-20^{\circ}$ providing positioning flexibility for extremity imaging and other challenging exams. The precise counterbalanced mechanics allows for an effortless and smooth travel of 56.5" from the floor. Soft touch, antimicrobial PA and multi-positional overhead handgrips aide in patient positioning and stability, while the flexible vertical drive positioning handle and conveniently located electromagnetic lock releases reduce technologist wrist strain and fatigue.

Key Highlights

Receptor Tilt Range: $+90^{\circ}/-20^{\circ}$

Available with Cabinet Options

- Standard cassette tray, rotating tray, or 17" x 17" fixed detector mounted inside

Smooth, Counterbalanced Movements

- Accommodating to all users regardless of size and strength

Electromagnetic Locks

- Conveniently placed control buttons are positioned on the side of the receptor cabinet

Vertical Travel Range of 56"

- Provides patient comfort for standing or seated chest examinations

Flexible Install

- Control buttons and load side can be changed to either side of the image receptor

Optional Accessories


- Overhead patient handgrips
- 8" extension box
- Rotating cassette tray

Technical Specifications

* All specifications subject to change without notice

VT300 Wall Stand

Receptor tilt range:	+90°/ 20°
Vertical travel range, vertical receptor:	15.5" (39.4 cm) to 71.1" (180.6 cm)
Vertical travel range, +90° receptor:	Minimum 25.5" (64.8 cm)
Height:	83.1" (211 cm)
Depth, vertical receptor:	30" (76 cm)
Receptor width:	24.6" (62.5 cm)
Receptor height:	26.5" (67.3 cm)
Weight:	765.5 lbs. (347.2 kg)
Weight, without counterweights:	368.5 lbs. (167 kg)
Power requirements:	120 VAC 230/240 VAC
Options:	- 8" (20 cm) extension box - Overhead patient handgrips - Rotating cassette tray

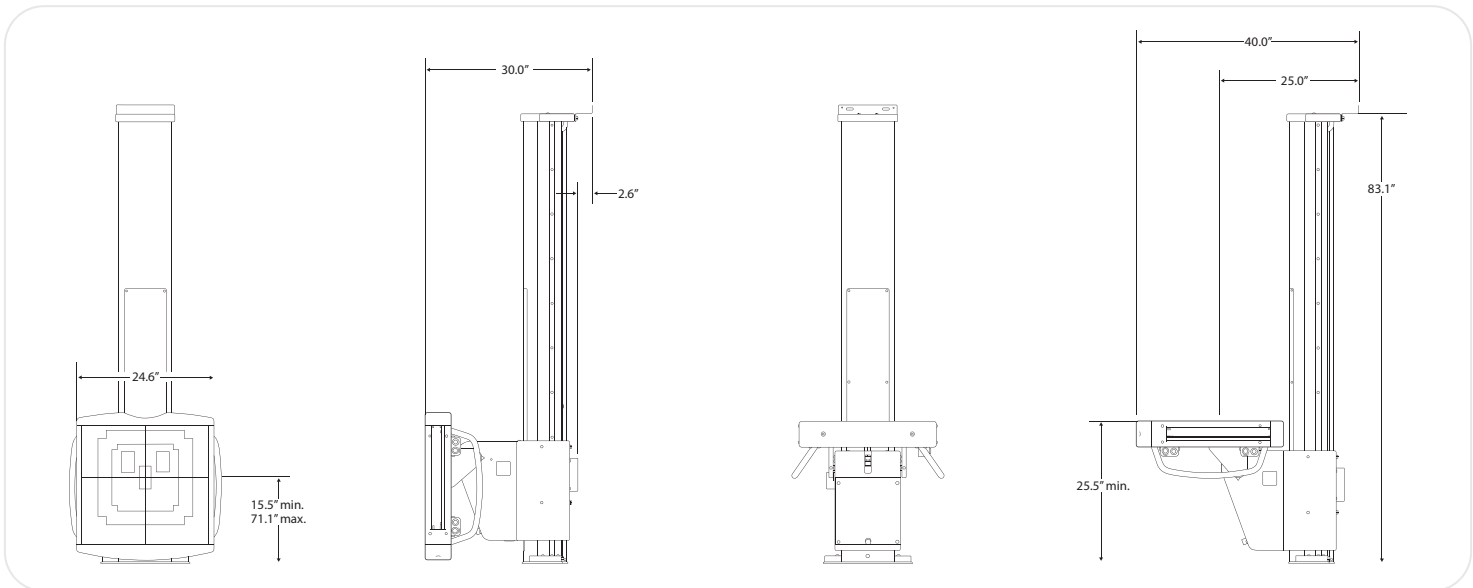

www.delmedical.com

New York
28 Calvert Street
Harrison, NY 10528
(Tel) 800.261.9808
(Fax) 914.835.6111

Illinois
241 Covington Drive
Bloomington, IL 60108
(Tel) 800.800.6006
(Fax) 847.288.7011

DEL MEDICAL

#9000-VT300 02.19